

Smart Guide KIT

The Smart Guide Kit provides a streamlined, efficient method for creating **customized surgical drill guides directly in the clinic**. The pre-formed templates are softened by a 70 °C water bath, molded to the patient's model using guide pins, and then used for precise implant drilling. It's a versatile, **system-independent** guide - ideal for standard clinical situations when conventional guided surgery isn't an option.



ADJUSTABLE & CHAIRSIDE FABRICATION

Thermoplastic templates soften in just 1 minute at 70 °C, allowing chairside customization for immediate surgical guidance.



MULTIPLE PRE-FORMED TEMPLATES

Available in four shapes (single bridge, short bridge, long bridge, and free-end bridge) to match common clinical indications.



SYSTEM-INDEPENDENT

Templates accommodate various implant systems, offering flexibility across different surgical setups.



Surgical Sequence



SmartGuide



Soak in warm water



Easy handling



Apply into the mouth



Drilling

Smart Guide KIT (OSGK)

Lower panel components:

4x Guide Pin
SGP22



2x Round bur
RAHM1018




2x Twist Drill (Ø2.2) For working models
2D2208LC01



SGCB30S

X Sleeve
(4ea)

Twist Drill (Ø2.2)
SGTD2207S

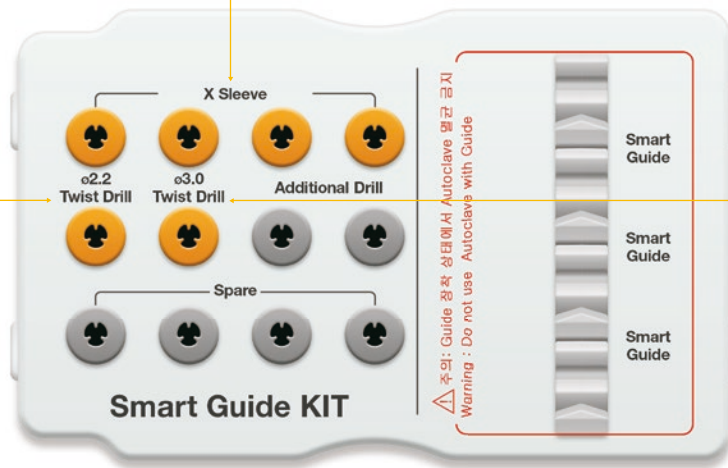

Twist Drill (Ø3.0)
QGTD3008

